


Aquaporin 1 Antibody / AQP1 (R32050)


Catalog No.	Formulation	Size
R32050	0.5mg/ml if reconstituted with 0.2ml sterile DI water	100 ug

Bulk quote request


Availability	1-3 business days
Species Reactivity	Human, Mouse, Rat
Format	Antigen affinity purified
Clonality	Polyclonal (rabbit origin)
Isotype	Rabbit IgG
Purity	Antigen affinity
Buffer	Lyophilized from 1X PBS with 2% Trehalose
UniProt	P29972
Localization	Membrane
Applications	Western Blot : 0.5-1ug/ml Immunohistochemistry (FFPE) : 2-5ug/ml Immunofluorescence (FFPE) : 5ug/ml Flow Cytometry : 1-3ug/million cells
Limitations	This Aquaporin 1 antibody is available for research use only.

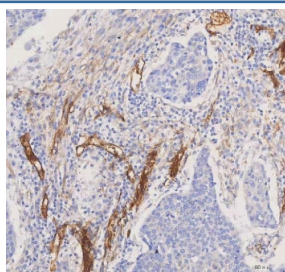


IHC staining of FFPE human bladder urothelial carcinoma tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

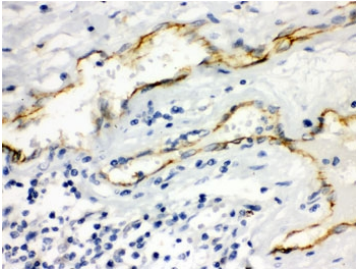


IHC testing of FFPE human intestinal cancer with Aquaporin 1 antibody. HIER: Boil the paraffin sections in pH 6, 10mM citrate buffer for 20 minutes and allow to cool prior to staining.

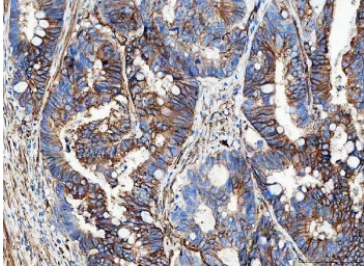


IHC staining of FFPE human rectal cancer tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

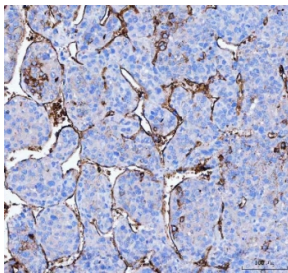


IHC staining of FFPE human liver cancer tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

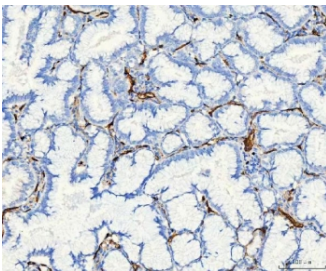


IHC staining of FFPE human lung cancer tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

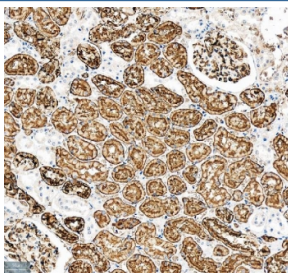


IHC staining of FFPE human kidney tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

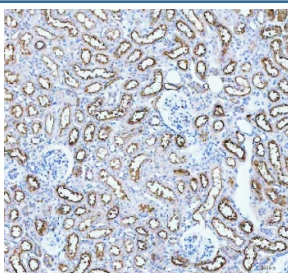


IHC staining of FFPE mouse kidney tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

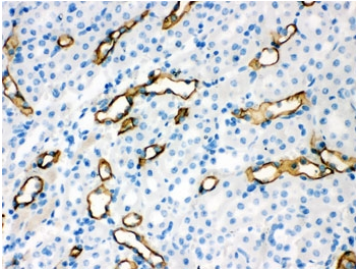

IHC testing of FFPE mouse kidney with Aquaporin 1 antibody. HIER: Boil the paraffin sections in pH 6, 10mM citrate buffer for 20 minutes and allow to cool prior to staining.

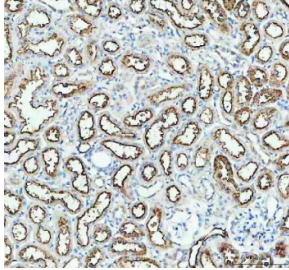

IHC staining of FFPE rat kidney tissue with Aquaporin 1 antibody. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

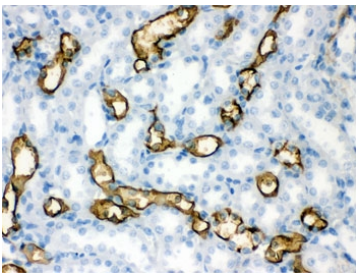

IHC testing of FFPE rat kidney with Aquaporin 1 antibody. HIER: Boil the paraffin sections in pH 6, 10mM citrate buffer for 20 minutes and allow to cool prior to staining.

Description

Aquaporin 1 is a 28-kD integral protein thought at first to be a breakdown product of the Rh polypeptide but was later shown to be a unique molecule that is abundant in erythrocytes and renal tubules. AQP1 is also expressed by the choroid plexus and various other tissues. It forms a water-specific channel that provides the plasma membranes of red cells and kidney proximal tubules with high permeability to water, thereby permitting water to move in the direction of an osmotic gradient.

Application Notes

Optimal dilution of the Aquaporin 1 antibody should be determined by the researcher.

Immunogen

Amino acids DRVKVWTSGQVEEYDLDDADDINSRVEMKPK of human Aquaporin 1 were used as the immunogen for the Aquaporin 1 antibody.

Storage

After reconstitution, the Aquaporin 1 antibody can be stored for up to one month at 4°C. For long-term, aliquot and store at -20°C. Avoid repeated freezing and thawing.

