


CD34 Antibody [clone SPM610] (V3001)


| Catalog No. | Formulation | Size |
|----------------|---|--------|
| V3001-100UG | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 100 ug |
| V3001-20UG | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 20 ug |
| V3001SAF-100UG | 1 mg/ml in 1X PBS; BSA free, sodium azide free | 100 ug |
| V3001IHC-7ML | Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only* | 7 ml |

Bulk quote request

| Availability | 1-3 business days |
|--------------------|--|
| Species Reactivity | Human |
| Format | Purified |
| Clonality | Monoclonal (mouse origin) |
| Isotype | Mouse IgG1, kappa |
| Clone Name | SPM610 |
| Purity | Protein G affinity chromatography |
| UniProt | P28906 |
| Localization | Cell surface |
| Applications | Immunohistochemistry (FFPE): 2-4ug/ml for 30 min at RT |
| Limitations | This CD34 antibody is available for research use only. |

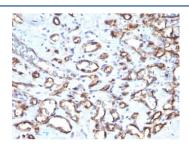

IHC: Formalin-fixed, paraffin-embedded human angiosarcoma stained with CD34 antibody (SPM610)

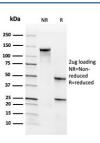

SDS-PAGE analysis of purified, BSA-free CD34 antibody (clone SPM610) as confirmation of integrity and purity.

Description

This antibody recognizes a carbohydrate epitope on a single chain, transmembrane, heavily glycosylated protein of 90-120kDa, which is identified as CD34 (VI international workshop on human differentiation antigens). Its expression is a hallmark for identifying pluripotent hematopoietic stem or progenitor cells. Its expression is gradually lost as lineage committed progenitors differentiate. CD34 is a marker of choice for staining blasts in acute myeloid leukemia. In addition, it is expressed by soft tissue tumors, such as solitary fibrous tumor and gastrointestinal stromal tumor. CD34 expression is also found in vascular endothelium. Additionally, proliferating endothelial cells overexpress this molecule than the non-proliferating endothelial cells. Anti-CD34 labels > 85% of angiosarcoma and Kaposi s sarcoma, but shows low specificity.

Application Notes

Optimal dilution of the CD34 antibody should be determined by the researcher.

- 1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 10-20 min followed by cooling at RT for 20 min.
- 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Immunogen

Recombinant full-length human protein was used as the immunogen for the CD34 antibody.

Storage

Store the CD34 antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).