


# CD45RB Antibody [clone PTPRC/1147] (V2415)

Catalog No.	Formulation	Size
V2415-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V2415-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V2415SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V2415IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml


### **Bulk quote request**


Species Reactivity	Human
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG1, kappa
Clone Name	PTPRC/1147
Purity	Protein G affinity chromatography
Buffer	1X PBS, pH 7.4
UniProt	P08575
Localization	Cell surface and cytoplasmic
Applications	Flow Cytometry: 0.5-1ug/10^6 cells Immunofluorescence: 0.5-1ug/ml Western Blot: 0.5-1ug/ml Immunohistochemistry (FFPE): 0.5-1ug/ml for 30 min at RT
Limitations	This CD45RB antibody is available for research use only.

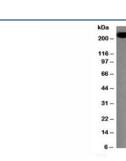

Western blot testing of CD45RB antibody and Daudi lysate. Expected molecular weight:  $\sim$ 147/180-220 kDa (unmodified/glycosylated).



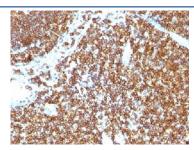

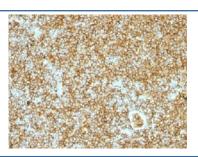

IHC testing of CD45RB antibody and FFPE human lymphoma

#### **Description**

CD45R, also designated CD45 and PTPRC, has been identified as a transmembrane glycoprotein, broadly expressed among hematopoietic cells. Multiple isoforms of CD45R are distributed throughout the immune system according to cell type. These isoforms arise because of alternative splicing of exons 4, 5, and 6. The corresponding protein domains are characterized by the binding of monoclonal antibodies specific for CD45RA (exon 4), CD45RB (exon 5), CD45RC (exon 6) and CD45RO (exons 4 to 6 spliced out). The variation in these isoforms is localized to the extracellular domain of CD45R, while the intracellular domain is conserved. CD45RB is expressed on mature B-lymphocytes and the majority of lymphomas and leukemias of B-cell origin. This antigen is also expressed on the surface of cells of myeloid lineage.

#### **Application Notes**

Due to variation in protocol and secondary antibody used, the CD45RB antibody may need to be titered for optimal performance.

- 1. FFPE staining requires boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes.
- 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

#### **Immunogen**

Human B lymphocytes were used as immunogen for this CD45RB antibody.

#### **Storage**

The CD45RB antibody (with azide) can be stored at 2-8oC. The azide-free format should be aliquoted and stored at -20oC or colder.

## References (2)