


Desmin Antibody [clone DES/3255] (V9390)


Catalog No.	Formulation	Size
V9390-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide	100 ug
V9390-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide	20 ug
V9390SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug

[Bulk quote request](#)


Availability	1-3 business days
Species Reactivity	Human, Dog, Cat, Horse, Cow
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG2b, kappa
Clone Name	DES/3255
Purity	Protein A/G affinity
UniProt	P17661
Localization	Cytoplasm
Applications	Immunohistochemistry (FFPE) : 2-4ug/ml
Limitations	This Desmin antibody is available for research use only.

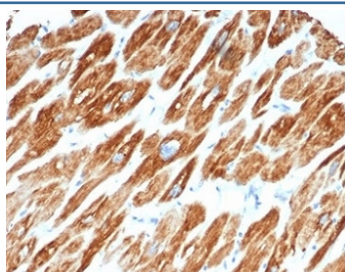


IHC staining of FFPE human cardiac muscle tissue with Desmin antibody (clone DES/3255) at 2ug/ml. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

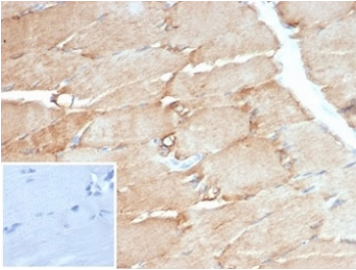


IHC staining of FFPE human skeletal muscle with Desmin antibody (clone DES/3255) at 2ug/ml. Negative control inset: PBS used instead of primary antibody to control for secondary Ab binding. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

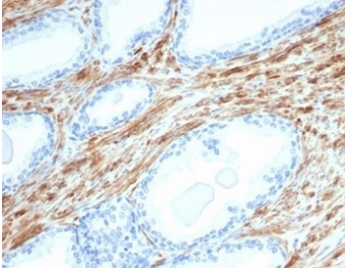


IHC staining of FFPE human prostate tissue with Desmin antibody (clone DES/3255) at 2ug/ml. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

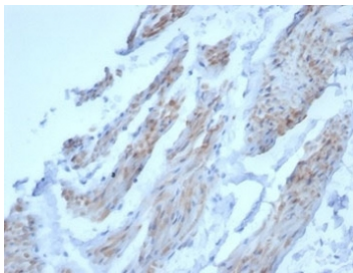


IHC staining of FFPE dog smooth muscle tissue with Desmin antibody (clone DES/3255) at 2ug/ml. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

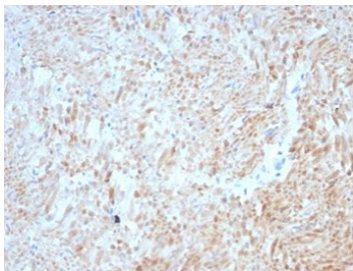


IHC staining of FFPE cat smooth muscle tissue with Desmin antibody (clone DES/3255) at 2ug/ml. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

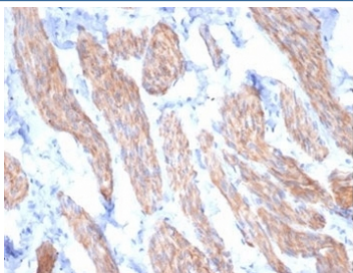

IHC staining of FFPE horse smooth muscle tissue with Desmin antibody (clone DES/3255) at 2ug/ml. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

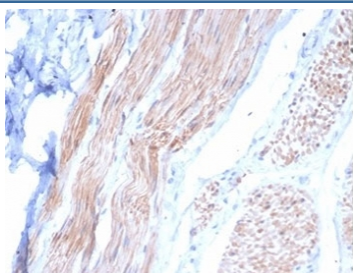

IHC staining of FFPE cow small intestine tissue with Desmin antibody (clone DES/3255) at 2ug/ml. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

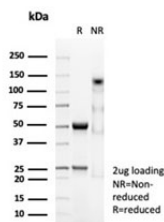

SDS-PAGE analysis of purified, BSA-free Desmin antibody (clone DES/3255) as confirmation of integrity and purity.

Description

Cytoskeletal intermediate filaments (IFs) constitute a diverse group of proteins that are expressed in a highly tissue-specific manner. IFs are constructed from two-chain α -helical coiled-coil molecules arranged on an imperfect helical lattice, and have been widely used as markers for distinguishing individual cell types within a tissue and identifying the origins of metastatic tumors. Vimentin is an IF general marker of cells originating in the mesenchyme. Vimentin and Desmin, a related class III IF, are both expressed during skeletal muscle development. Desmin, a 469 amino acid protein found near the Z line in sarcomeres, is expressed more frequently in adult differentiated state tissues. Anti-desmin detects cells of normal smooth, skeletal, and cardiac muscles. Antibody reacts with leiomyomas, leiomyosarcoma, rhabdomyomas, rhabdomyosarcoma, and perivascular cells of glomus tumors of the skin.

Application Notes

Optimal dilution of the Desmin antibody should be determined by the researcher.

Immunogen

Recombinant full-length human desmin protein was used as the immunogen for the Desmin antibody.

Storage

Aliquot the Desmin antibody and store frozen at -20°C or colder. Avoid repeated freeze-thaw cycles.