


hCG Holo Antibody [clone HCGab/52] (V3064)


Catalog No.	Formulation	Size
V3064-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V3064-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V3064SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V3064IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml

Bulk quote request

Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG1, kappa
Clone Name	HCGab/52
Purity	Protein G affinity chromatography
UniProt	P01215, P0DN86
Localization	Cytoplasmic, secreted
Applications	Immunohistochemistry (FFPE) : 0.5-1ug/ml for 30 min at RT
Limitations	This hCG Holo antibody is available for research use only.

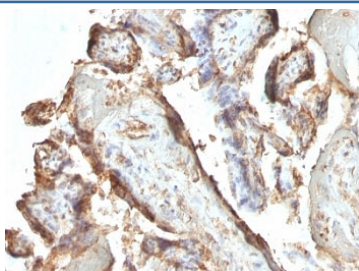

IHC: Formalin-fixed, paraffin-embedded human placenta stained with hCG Holo antibody (clone HCGab/52).

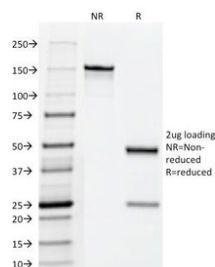

SDS-PAGE Analysis of Purified, BSA-Free hCG Holo Antibody (clone HCGab/52).
Confirmation of Integrity and Purity of the Antibody.

Description

This mAb reacts only with the intact-HCG and not with either free alpha- or free beta-chain of HCG. HCG is a glycoprotein and is composed of two non-identical, non-covalently linked polypeptide chains designated as the alpha and beta subunits. The alpha subunit is identical to that of thyroid stimulating hormone (TSH), follicle stimulating hormone (FSH), and luteinizing hormone (LH). HCG is secreted in large quantities by normal trophoblasts. It is present only in trace amounts in non-pregnant urine and sera but rises sharply during pregnancy. HCG mAb detects cells and tumors of trophoblastic origin such as choriocarcinoma. Large cell carcinoma and adenocarcinoma of the lung demonstrate anti-hCG positivity in 90% and 60% of cases respectively. 20% of lung squamous cell carcinomas are positive. HCG expression by non-trophoblastic tumors may indicate aggressive behavior.

Application Notes

Optimal dilution of the hCG Holo antibody should be determined by the researcher.

1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min.
2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Immunogen

Purified hCG protein was used as the immunogen for the hCG Holo antibody.

Storage

Store the hCG Holo antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).