


PLAP Antibody / Placental Alkaline Phosphatase [clone SPM593] (V2515)

Catalog No.	Formulation	Size
V2515-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V2515-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V2515SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V2515IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml

Bulk quote request

Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Clonality	Monoclonal (mouse origin)
Isotype	Mouse IgG2b, kappa
Clone Name	SPM593
Purity	Protein G affinity chromatography
UniProt	P05187
Localization	Cytoplasmic and cell surface
Applications	Immunohistochemistry (FFPE) : 1-2ug/ml for 30 min at RT
Limitations	This PLAP antibody is available for research use only.

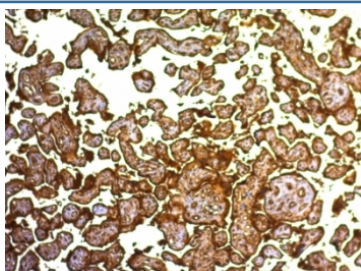

IHC: Formalin-fixed, paraffin-embedded human placenta stained with PLAP antibody (clone SPM593).

Description

Reacts with a 70kDa membrane-bound isozyme (Regan and Nagao type) of placental Alkaline Phosphatase (PLAP) occurring in the placenta during the 3rd trimester of gestation. It is highly specific for PLAP and shows no cross-reaction with other isozymes of alkaline phosphatase. Anti-PLAP reacts with germ cell tumors and can discriminate between these and other neoplasms. Somatic neoplasms e.g. breast, gastrointestinal, prostatic, and urinary cancers may also immunoreact with antibodies to PLAP. Anti-PLAP positivity in conjunction with anti-keratin negativity favors seminoma over carcinoma. Germ cell tumors are usually anti-keratin positive, but they regularly fail to stain with anti-EMA, whereas most carcinomas stain with anti-EMA. Anti-PLAP has been useful in the diagnosis of gestational trophoblastic disease.

Application Notes

Optimal dilution of the PLAP antibody should be determined by the researcher.

1. No special pretreatment is required for the immunohistochemical staining of formalin-fixed tissues.
2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Immunogen

Recombinant full-length human ALPP/placental alkaline phosphatase protein was used as the immunogen for the PLAP antibody.

Storage

Store the PLAP antibody at 2-8°C (with azide) or aliquot and store at -20°C or colder (without azide).