


# Recombiant Ki-67 Antibody [clone MKI67/4945R] (V9265)

Catalog No.	Formulation	Size
V9265-100UG	0.2~mg/ml in 1X PBS with $0.1~mg/ml$ BSA (US sourced), $0.05%$ sodium azide	100 ug
V9265-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced), 0.05% sodium azide	20 ug
V9265SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug


### Recombinant RABBIT MONOCLONAL

### **Bulk quote request**


Availability	1-3 business days
Species Reactivity	Human
Format	Purified
Clonality	Recombinant Rabbit Monoclonal
Isotype	Rabbit IgG
Clone Name	MKI67/4945R
Purity	Protein A/G affinity
UniProt	P46013
Localization	Nuclear
Applications	Western Blot : 1-2ug/ml Immunohistochemistry (FFPE) : 1-2ug/ml
Limitations	This recombinant Ki-67 antibody is available for research use only.

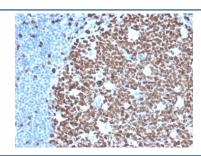


IHC staining of FFPE human tonsil tissue with recombinant Ki-67 antibody (clone MKI67/4945R). HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

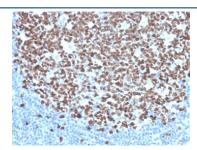


IHC staining of FFPE human tonsil tissue with recombinant Ki-67 antibody (clone MKI67/4945R). HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

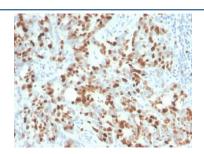

IHC staining of FFPE human lymph node tissue with recombinant Ki-67 antibody (clone MKI67/4945R). HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

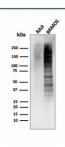

Western blot testing of human Raji and Ramos cell lysates using recombinant Ki-67 antibody (clone MKI67/4945R).

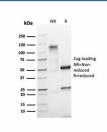

SDS-PAGE analysis of purified, BSA-free recombinant Ki-67 antibody (clone MKI67/4945R) as confirmation of integrity and purity.

## **Description**

Ki-67 antigen is a nuclear, non-histone protein that is present in all stages of the cell cycle except G0. This characteristic makes Ki-67 an excellent marker for proliferating cells and is commonly used as one of the prognostic factors in cancer studies. A correlation has been demonstrated between Ki-67 index and the histo-pathological grade of neoplasms. Assessment of Ki-67 expression in renal and ureter tumors shows a correlation between tumor proliferation and disease progression, thus making it possible to differentiate high-risk patients. Ki-67 expression may also prove to be important for distinguishing between malignant and benign peripheral nerve sheath tumors. Ki-67 labeling index has been shown to be a prognostic marker in a number of neoplasms including grade II astrocytoma, oligodendroglioma, colon carcinoma, and breast carcinoma. In general, Ki-67 is a good marker of proliferating cell populations.

# **Application Notes**

Optimal dilution of the recombinant Ki-67 antibody should be determined by the researcher.

### **Immunogen**

A portion of amino acids 2293-2478 from the human protein was used as the immunogen for the recombinant Ki-67 antibody.

# **Storage** Aliquot the recombinant Ki-67 antibody and store frozen at -20oC or colder. Avoid repeated freeze-thaw cycles.