


# Recombinant IgM Heavy Chain Antibody [clone rlM373] (V7425)

Catalog No.	Formulation	Size
V7425-100UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	100 ug
V7425-20UG	0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide	20 ug
V7425SAF-100UG	1 mg/ml in 1X PBS; BSA free, sodium azide free	100 ug
V7425IHC-7ML	Prediluted in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide; *For IHC use only*	7 ml


# Recombinant MOUSE MONOCLONAL

### **Bulk quote request**

Species Reactivity	Human
Format	Purified
Clonality	Recombinant Mouse Monoclonal
Isotype	Mouse IgG1, kappa
Clone Name	rlM373
Purity	Protein G affinity chromatography
Gene ID	3507
Localization	Cytoplasm, cell surface and secreted
Applications	Immunohistochemistry (FFPE): 0.5-1ug/ml for 30 min at RT
Limitations	This recombinant IgM Heavy Chain antibody is available for research use only.

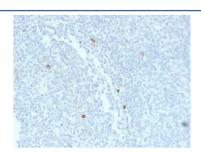

IHC staining of FFPE human tonsil tissue with recombinant IgM Heavy Chain antibody (clone rIM373). HIER: boil tissue sections in pH6, 10mM citrate buffer, for 10-20 min and allow to cool before testing.

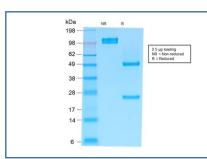

SDS-PAGE analysis of purified, BSA-free recombinant IgM Heavy Chain antibody (clone rlM373) as confirmation of integrity and purity.

## **Description**

This antibody recognizes a protein of 75kDa, identified as the mu heavy chain of human immunoglobulins. It does not cross-react with alpha (IgA), gamma (IgG), epsilon (IgE), or delta (IgD), heavy chains, T-cells, monocytes, granulocytes, or erythrocytes. The antibody is useful in the identification of leukemias, plasmacytomas, and certain non-Hodgkin's lymphomas. The most common feature of these malignancies is the restricted expression of a single heavy chain class. Demonstration of clonality in lymphoid infiltrates indicates that the infiltrate is clonal and therefore malignant.

## **Application Notes**

The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the recombinant Heavy Chain antibody to be titered up or down for optimal performance.

1. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

#### **Immunogen**

Recombinant protein was used as the immunogen for this recombinant IgM Heavy Chain antibody.

#### **Storage**

Store the recombinant IgM Heavy Chain antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).

References (1)