


# TPO Antibody / Thyroid Peroxidase [clone TPO/3701] (V8625)


| Catalog No.    | Formulation  | Size   |
|----------------|--|--------|
| V8625-100UG    | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 100 ug |
| V8625-20UG     | 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide | 20 ug  |
| V8625SAF-100UG | 1 mg/ml in 1X PBS; BSA free, sodium azide free                             | 100 ug |

## **Bulk quote request**


| Availability       | 1-3 business days  |
|--------------------|--|
| Species Reactivity | Human  |
| Format             | Purified   |
| Clonality          | Monoclonal (mouse origin)                                  |
| Isotype            | Mouse IgG2b, kappa   |
| Clone Name         | TPO/3701   |
| Purity             | Protein G affinity chromatography                          |
| UniProt            | P07202   |
| Localization       | Cytoplasmic, cell surface                                  |
| Applications       | Immunohistochemistry (FFPE): 1-2ug/ml for 30 minutes at RT |
| Limitations        | This TPO antibody is available for research use only.      |

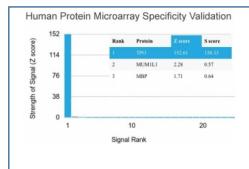

Analysis of HuProt(TM) microarray containing more than 19,000 full-length human proteins using TPO antibody. These results demonstrate the foremost specificity of the TPO/3701 mAb. Z- and S- score: The Z-score represents the strength of a signal that an antibody (in combination with a fluorescently-tagged anti-IgG secondary Ab) produces when binding to a particular protein on the HuProt(TM) array. Z-scores are described in units of standard deviations (SD's) above the mean value of all signals generated on that array. If the targets on the HuProt(TM) are arranged in descending order of the Z-score, the S-score is the difference (also in units of SD's) between the Z-scores. The S-score therefore represents the relative target specificity of an Ab to its intended target.

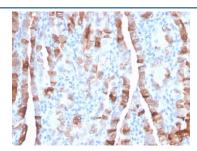

IHC staining of FFPE human thyroid carcinoma with TPO antibody. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

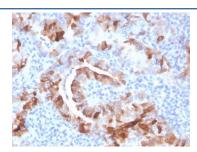

IHC staining of FFPE human thyroid carcinoma with TPO antibody. HIER: boil tissue sections in pH 9 10mM Tris with 1mM EDTA for 20 min and allow to cool before testing.

# **Description**

Thyroperoxidase (TPO) is a 933 amino acid, type I transmembrane glycoprotein that plays a key role in thyroid gland function and autoimmunity. It is present as a dimer on the apical surface of thyroid follicular cells. TPO functions in the iodination of tyrosine residues in thyroglobulin and phenoxy-ester formation between pairs of iodinated tyrosines to generate the thyroid hormones, thyroxine and triiodothyronine. Mutations in this gene are associated with several disorders of thyroid hormonogenesis, including congenital hypothyroidism, congenital goiter, and thyroid hormone organification defect IIA.Malignant thyroid tumors exhibit an anomaly in TPO resulting in lower affinity for anti-TPO. This antibody may aid in the differentiation between benign and malignant thyroid tumors.

### **Application Notes**

Optimal dilution of the TPO antibody should be determined by the researcher.

#### **Immunogen**

A portion of amino acids 685-804 from the human protein was used as the immunogen for the TPO antibody.

#### **Storage**

Store the TPO antibody at 2-8oC (with azide) or aliquot and store at -20oC or colder (without azide).