


Zebrafish Leo1 Antibody / RNA polymerase-associated protein LEO1 (RZ1233)


Catalog No.	Formulation	Size
RZ1233	0.5mg/ml if reconstituted with 0.2ml sterile DI water	100 ug

[Bulk quote request](#)


Availability	2-3 weeks
Species Reactivity	Zebrafish
Format	Antigen affinity purified
Clonality	Polyclonal (rabbit origin)
Isotype	Rabbit Ig
Purity	Antigen affinity chromatography
Buffer	Lyophilized from 1X PBS with 2% Trehalose
UniProt	Q6NYV9
Localization	Nuclear
Applications	Western Blot : 0.5-1ug/ml Immunohistochemistry (FFPE) : 2-5ug/ml Immunofluorescence : 5ug/ml
Limitations	This Zebrafish Leo1 antibody is available for research use only.

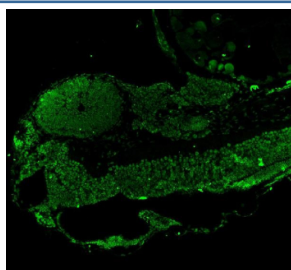


Immunofluorescent staining of Leo1 protein using Zebrafish Leo1 antibody (green) and zebrafish embryo tissue. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

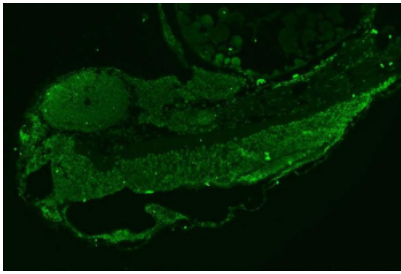


Immunofluorescent staining of Leo1 protein using Zebrafish Leo1 antibody (green) and zebrafish embryo tissue. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

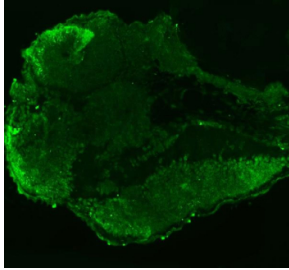

Immunofluorescent staining of Leo1 protein using Zebrafish Leo1 antibody (green) and zebrafish embryo tissue. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

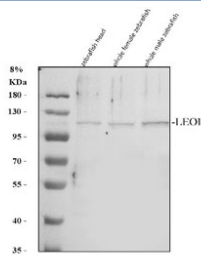


Western blot analysis of Leo1 protein using Zebrafish Leo1 antibody and 1) zebrafish head, 2) whole female zebrafish and 3) whole male zebrafish tissue lysate. Predicted molecular weight ~76 kDa but can be observed at ~105 kDa (human similarity).

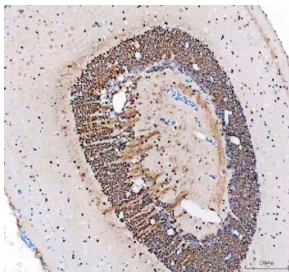

IHC staining of FFPE zebrafish brain tissue with Zebrafish Leo1 antibody, HRP-labeled secondary and DAB substrate. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.



IHC staining of FFPE zebrafish eye tissue with Zebrafish Leo1 antibody, HRP-labeled secondary and DAB substrate. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

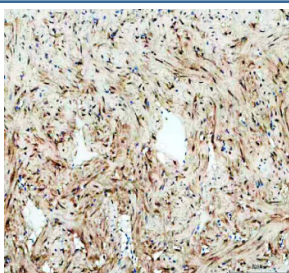

IHC staining of FFPE zebrafish heart tissue with Zebrafish Leo1 antibody, HRP-labeled secondary and DAB substrate. HIER: boil tissue sections in pH8 EDTA for 20 min and allow to cool before testing.

